

Using a Small Volume Compressor Nebulizer



1. Put the equipment on a level, steady surface. Plug the power cord into a wall outlet.



2. Wash your hands. Connect one end of the tube to the air outlet connector.



3. Add the medicine as prescribed with the unit dose vial.




4. Connect the mouthpiece, T-piece and reservoir tubing to the nebulizer.



5. Attach the tubing to the air inlet.



6. Turn on the power switch.

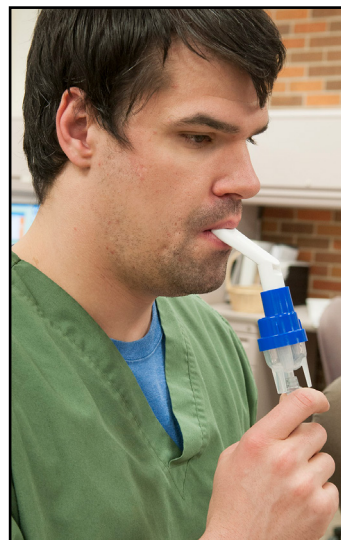

7.

Sit in a comfortable position. Put the mouthpiece between your upper and lower teeth, over your tongue and with your lips closed.

Breathe normally, but take a deeper breath every minute until the nebulizer begins to sputter.

Turn off the power switch.

(over)

Cleaning a Small Volume Compressor Nebulizer



1. Disconnect the nebulizer from the tubing. Take all the pieces apart. Wash in warm, soapy water.



2. Rinse the nebulizer parts in warm water and gently shake off any water.



3. Lay the parts on a clean towel or dish rack to air dry. Once dry, put the pieces together and keep in a cool, dry place.